


## Movettes Pipettes

- Made of polyethylene
- 3 sizes : 1.5, 4 and 5 ml
- For semi-quantitative sampling

Cat. No.	Volume (ml)	Sterile	Units/carton	€ Excl. VAT/carton
391474*	1.5	No	100	NC -
391475*	5	No	100	NC -

\* Economic models

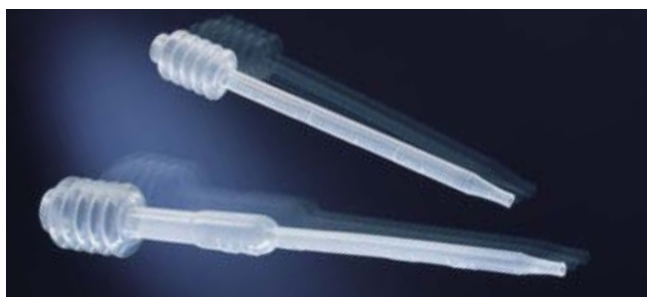


## Wide-mouth transfer pipettes

Facilitate the handling of fragile or high viscosity products.

- In LDPE
- Large bulb
- Total length : 125 mm
- Internal diameter of the tip : 5.5 mm

Cat. No.	Description	Units/box	€ Excl. VAT/box
380514	Wide-mouth transfer pipette	500	NC -

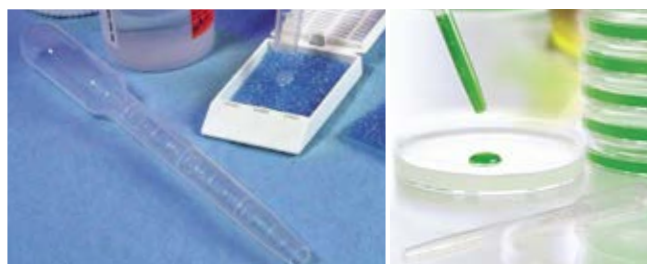


## Fast-Read 102® Plastic counting cells

10-grid PPMA slides for counting cells per  $\mu$ l.

Cell volume ( $\mu$ l)	7
Grid dimensions (mm)	5 x 2
Depth (mm)	0.1
Grid volume ( $\mu$ l)	1
Square dimensions (mm)	1 x 1
Square volume ( $\mu$ l)	0.1
Chamber volume ( $\mu$ l)	37

Cat. No.	Description	€ Excl. VAT
390497	Fast-read® slides, 10 chambers x 100	NC -
390498	Fast-Read® 100 slides Kit, 10 chambers + 1000 transfer pipettes	NC -
390499	Fast-Read® 100 slides kit, 10 chambers + 1000 transfer pipettes + 1 bottle of dye 35 ml	NC -

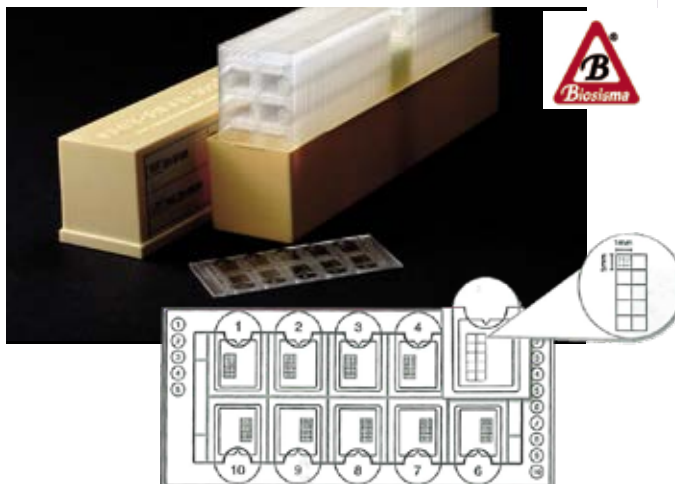

## M-T® System for Fast-Read 102®

The use of the M-T® system combined with the Fast-Read 102® slide considerably reduces the risk of confusion with other cells, thus increasing the accuracy of the result.

- Used for urine examinations, avoids overlapping of epithelial cells during sedimentation
- M-T® transfer pipettes: special calibration allowing 1 ml of sample to be blocked at the bottom of the tube
- Tubes for M-T® system: 16 x 103 mm, PS transparent, graduation 1/2,5/5/10 ml

Cat. No.	Description	€ Excl. VAT
390508	1000 M-T® transfer pipettes + 1000 x 10 ml tubes	NC -
391536	Rack For 10 tubes	NC -
390503	Kit of 100 Fast-read® 10-chamber slides + 1000 M-T® transfer pipettes	NC -
390506	Kit of 100 Fast-read® 10-chamber slides + 1000 M-T® transfer pipettes + 1000 x 10 ml tubes + 1000 caps	NC -
390522	Conical tubes 16 x 103 mm graduated in PS (2000 units)	NC -
390523	Conical tubes 16 x 103 mm graduated in PS (2000 units)	NC -
390517	Tri-bulb pipettes for M-T® system (3000 units)	NC -



### ➤ Principle

1. Place tube 390522 or 390523 containing the urine to be analysed in the express rack 391536 and centrifuge
2. Insert the 3-bulb M-T transfer pipette
3. Turn the rack over
4. You get 1 ml of resuspended sediment ready for analysis on slide Fast-Read 102®

